


Quick De-Gas

BQSL

Quick De-gas is an instrument that degasses a sample in an incredibly short time. Ingredients in a sample are influenced by time of degassing procedure, so it should be done in a short time. Quick De-gas completes this procedure efficiently in a short time with a combination of stirring, depressing, and aeration. For example, it takes only 3 minute to degas a sample of 200ml.

If you use saturation bottle, evaporation of high volatile ingredients such as alcohol can be also controlled.

This instrument enables you an efficient testing procedure by saving time.

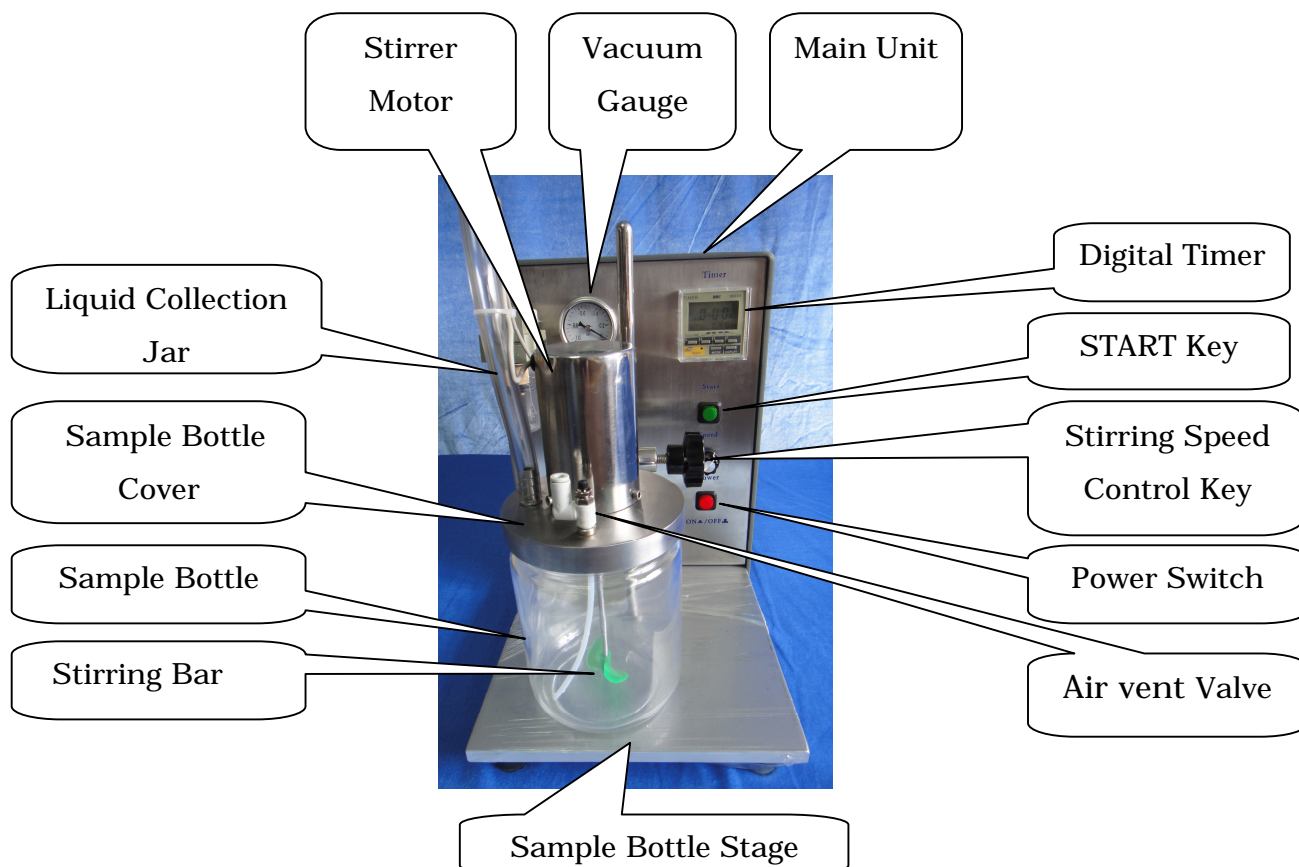

Specifications / Quick De-Gas

	Suction System
Model	QD11
Voltage (V)	AC 100-240
Frequency (Hz)	50/60
Power (mA)	2000
Vacuum (kPa)	-80
Capacity	38 liter/minute
Length (mm)	480
Width (mm)	180
Height (mm)	450
Weight (kg)	8.0
Time Set	0-99 minutes

Specifications and outer appearance are subject to change without notice.

Make your process simple!

Higher quality is expected!

For more information, please contact :

BQSL Japan Co., Ltd.

5-11-8-102 Okada, Atsugi-shi, Kanagawa,
243-0021 JAPAN

Tel / Fax: +81-46-230-0417

e-mail: mail@bqsl.com

<http://www.bqsl.com>